



# Ei2300 Monitoring System

## Siloxanes

Complex and varying gas matrices make siloxanes difficult and expensive to measure. The Ei2300 Siloxane Monitoring System is an independent and fully automated measurement system developed to precisely quantify the concentration of siloxanes in biogas. The instrument separates compounds in two separate stages, combining the selectivity of a Gas Chromatograph (GC) with the extraordinary sensitivity (low ppbv) of an Ion Mobility Spectrometer (IMS), and is now supported by ASTM Standard D8455.

### Applications

- Biogas
- Dairy, Hog and Pig Farms
- Gaseous Fuels
- Landfills
- Wastewater Treatment Plants
- Food Waste Digesters (co-digestion)
- Others

### Available Configurations

- Tabletop
- Rack-Mounted
- Commercial and Industrial Enclosure
  - NEMA 4/IP66 Rated
  - Table and Wall Mountable
  - Window to View LCD Screen
  - Full Remote Access
  - Secure Enclosure
  - Integrated, Heated Sampling System
  - Heated Sample Line
  - Options
    - Hazardous Rated Enclosure (Class 1, Division 2)
    - Side Mounted Air Conditioning Unit
    - Remote Access Hardware



Supported by Standard D8455



Supported by Standard  
2613-2:2023






## Key Features

- Precise Measurements
- Easy to Operate
- No Analytical Specialist Needed
- Two-fold Matrix Separation
  - Gas Chromatograph
  - Ion Mobility Spectrometer

Ohio Lumex Proprietary Calibration Procedure Typical Limits of Quantitation (LOQ) [1 mL Loop] *		
Species	LOQ (mg Si/m <sup>3</sup> )	LOQ (ppb)
(TMSOL) Trimethylsilanol	0.005	1.2
(L2) Hexamethyldisiloxane	0.005	0.8
(L3) Octamethyltrisiloxane	0.006	0.6
(L4) Decamethyltetrasiloxane	0.006	0.5
(L5) Dodecamethylpentasiloxane	0.006	0.4
(D3) Hexamethylcyclotrisiloxane	0.006	0.7
(D4) Octamethylcyclotetrasiloxane	0.007	0.6
(D5) Decamethylcyclopentasiloxane	0.007	0.5
(D6) Dodecamethylcyclohexasiloxane	0.008	0.4

\* Example, may vary by application

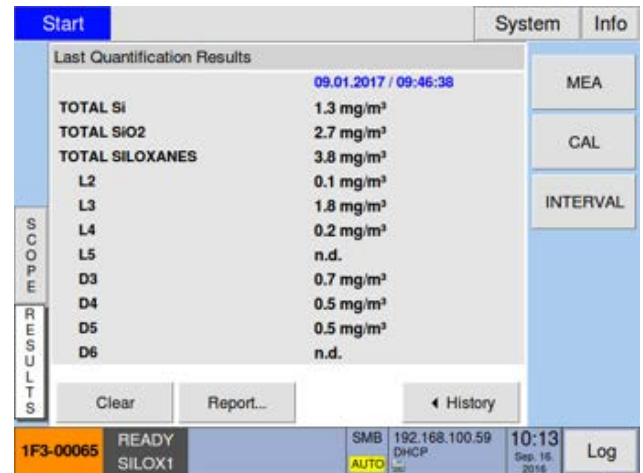

Last Quantification Results		System	Info
		09.01.2017 / 09:46:38	
TOTAL Si	1.3 mg/m <sup>3</sup>	MEA	
TOTAL SiO <sub>2</sub>	2.7 mg/m <sup>3</sup>	CAL	
TOTAL SILOXANES	3.8 mg/m <sup>3</sup>	INTERVAL	
L2	0.1 mg/m <sup>3</sup>		
L3	1.8 mg/m <sup>3</sup>		
L4	0.2 mg/m <sup>3</sup>		
L5	n.d.		
D3	0.7 mg/m <sup>3</sup>		
D4	0.5 mg/m <sup>3</sup>		
D5	0.5 mg/m <sup>3</sup>		
D6	n.d.		

Technical Specifications			
Separation I	GC Retention Time	Purge Gas Requirement	Nitrogen: Quantity 5.0 (cylinder or generator)
Gas Chromatograph (GC)	Isotherm (< 80° C)	Operation	Manual: 6.4" TFT touchscreen Online Remote Test: User defined intervals
Capillary Column	30m	Data Storage (Internal)	16 GB Flash Memory
Separation II	Drift Time Ion Mobility Spectrometer (TOF-IMS)	Data Output	USB, 4-20mA, MODBUS (TCP)
Ionization	β-radiation - <sup>3</sup> H	Ambient Temperature	Operation 0 - 40° C
Activity	300 MBq, below the exemption limit of 1 GBq acc. To EURATOM guideline	Moisture	Up to 95% non-condensing
Detection	Electrometer, Ion Mobility Spectrometer	Power Range	100 - 240 VAC, 50 - 60 Hz
LoD for Siloxanes	One digit ppb	Power Consumption	< 200 Watt
Typical Range for Siloxanes	0 - 1.5 mg/m <sup>3</sup>	Weight	34 lbs (15.5 kg)
Pneumatics	Two Electronic Pressure Controllers (EPC) for flow stability and GC ramping	Dimensions	17.7 x 17.1 x 7 in (No Case) 449 x 435 x 177 mm (No Case)
Sampling	Ambient pressure, Firmware controlled pump plus heated electrical six port valve, Liquid trap (water condensation) optional	Housing	19", IP 20 Enclosure, CE Marking

GC-IMS Supplied By:

